



**FIXART**  
SPORTS MEDICINE



# KNEE PORTFOLIO

# ACL / PCL Reconstruction

## Grapple Unity

Adjustable Loop Cortical Fixation Implant

### FEATURES

- Triple Loop for better graft fixation
- Loop Sleeve for graft preservation
- #5 UHMWPE Suture for maximum strength and smooth adjustment
- Triple assurance with mechanical loop lock, knot clutching and serial loop reduction
- Graft tensions achieved upto 900N
- Docking dome at bottom for secure positioning over the tunnel & arrest migration



## LigaZip B

Bioabsorbable Interference Screws

### FEATURES

- Dual Start thread for quick insertion
- Blunt buttress threads for maximum clinch without graft wear
- Made from  $\beta$ -TCP + PLA/PLGA mixture for optimised bone conductivity and strength




# LigaZip T

Titanium Interference Screws

## FEATURES

- Dual Start thread for quick insertion
- Blunt buttress threads for maximum clinch without graft wear

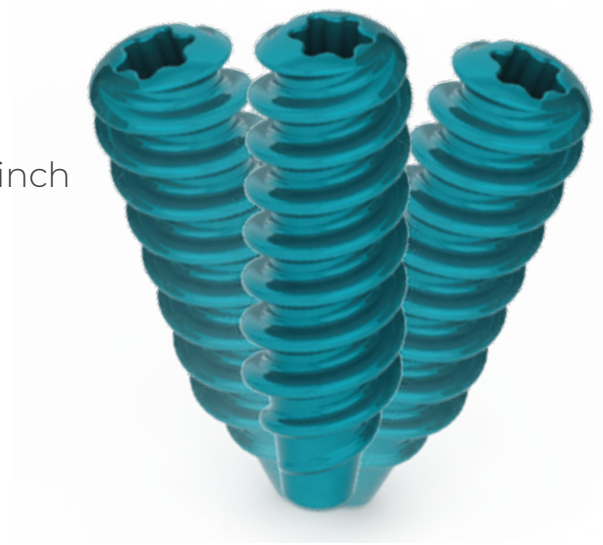

# LigaZip P

PEEK Interference screws

## FEATURES

- Dual Start thread for quick insertion
- Blunt buttress threads for maximum clinch without graft wear
- Made from PEEK-Optima for superior strength and Biocompatibility




# Meniscus Repair

## Menizip

All-Inside Meniscal Repair System

### FEATURES

- **Versatile All-Inside Repair:** Available in Dual loaded (1 Pair) implants and Quad loaded (2 Pair) implants variants
- **Effortless Single-Hand Deployment:** Ease of orientation and rotation as per comfort and approach, ensuring precise usability
- Superior strength and biocompatibility with PEEK Optima implants
- Achieve minimal incisions with a needle diameter of just 1.5 mm
- **Real-time Feedback:** Track the actual status of deployment by position marker for precision and flawless deploying

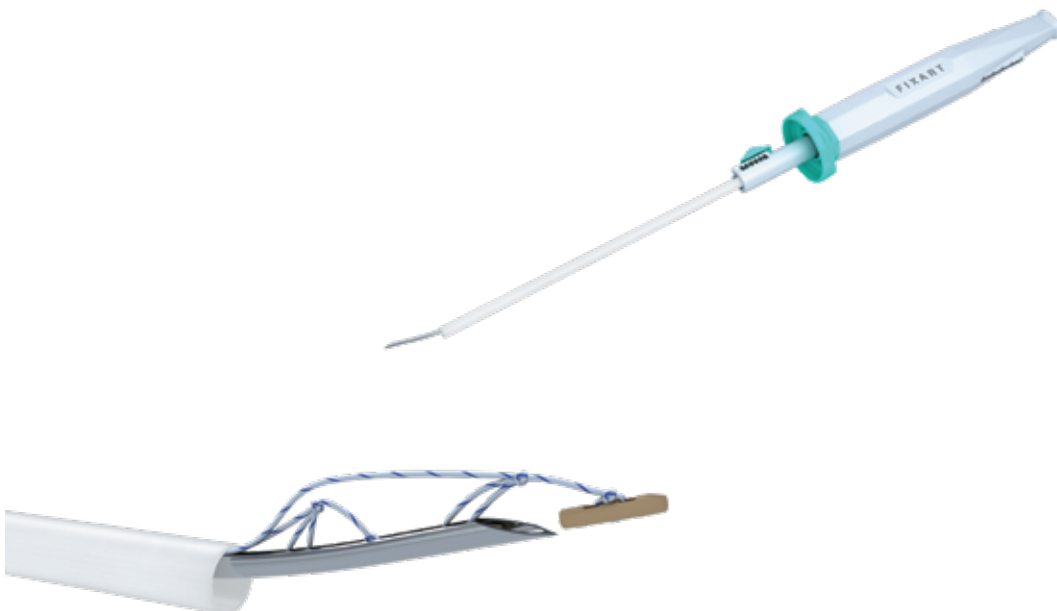



# SHOULDER PORTFOLIO

# Rotator Cuff Repair

## Deepfix

Screw-In PEEK Suture Anchor

### FEATURES

- Buttress thread design with dual start at cortical region for maximum pullout strength
- Loaded with #2 UHMWPE striped and coloured sutures
- Available in Double Loaded and Triple Loaded Options
- Made from PEEK-Optima for superior strength and biocompatibility



## CortiZip

Soft Tissue Repair PEEK Anchor

### FEATURES

- Eyelet designed for passing upto 6 #2 sutures
- Stationary tension retainment mechanism to maintain steady tension till full deployment
- Unique self holding distal push anchor at partial deployment with barbed design
- Mutual docking of split anchors at deployment is achieved for stability & alignment
- Made with PEEK-Optima for superior strength and biocompatibility







# Tenemorf RC

All-Suture Anchor

## FEATURES

- Self-Punching deployer
- Loaded with #2 UHMWPE Sutures
- Available in Double Loaded and Triple Loaded Options
- UHMWPE + PET fold soft anchor with subcortical deployment for maximum pullout strength

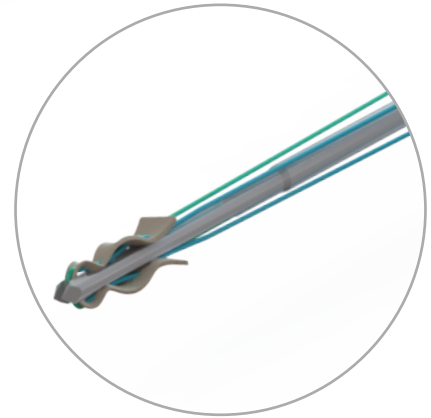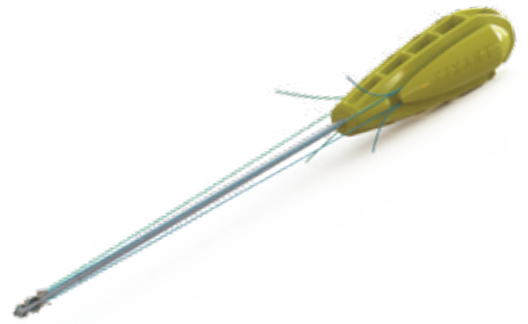

## Bankart/ Slap Repair

# PeriZip

PEEK Suture Anchor

## FEATURES

- Made with PEEK-Optima for superior strength and biocompatibility
- Loaded with #2 UHMWPE sutures
- Available in Single Loaded and Double Loaded Options





# Tenemorf

All-Suture Anchor

## FEATURES

- Flexible Deployer for Curved Guides and Sleeves
- Loaded with #2 UHMWPE Sutures
- Available in Single Loaded and Double Loaded Options
- UHMWPE + PET fold soft anchor with sub cortical deployment for maximum pullout strength

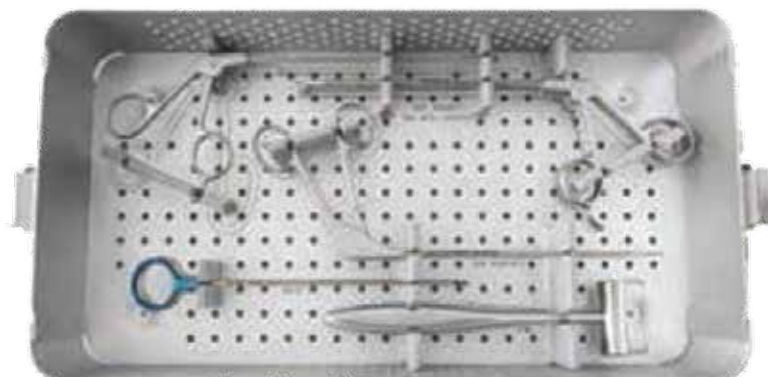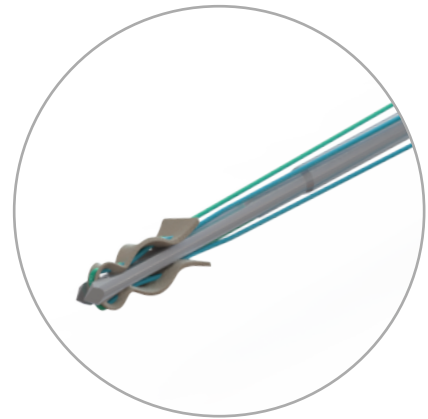

# Disposable Cannulas

Disposable Full Thread  
Arthroscopy Cannula



# Cat Numbers

## Grapple Unity

Adjustable Loop Cortical Fixation Implant

CAT NO.

3571

## LigaZip B

Bioabsorbable Interference Screws

CAT NO.	SIZES
3582	Ø7 x 25 mm
3583	Ø7 x 30 mm
3584	Ø8 x 25 mm
3585	Ø8 x 30 mm
3586	Ø9 x 25 mm
3587	Ø9 x 30 mm
3588	Ø10x 25 mm
3589	Ø10x 30 mm
3590	Ø11x 25 mm
3591	Ø11x 30 mm

## LigaZip T

Titanium Interference Screws

CAT NO.	SIZES
3572	Ø7 x 25 mm
3573	Ø7 x 30 mm
3574	Ø8 x 25 mm
3575	Ø8 x 30 mm
3576	Ø9 x 25 mm
3577	Ø9 x 30 mm
3578	Ø10x 25 mm
3579	Ø10x 30 mm
3580	Ø11x 25 mm
3581	Ø11x 30 mm

# Cat Numbers

## LigaZip P

PEEK Interference screws

CAT NO.	SIZES
7 x 25	4036
7 x 30	4037
8 x 25	4038
8 x 30	4039
9 x 25	4040
9 x 30	4041
10 x 25	4042
10 x 30	4043
11 x 25	4044
11 x 30	4045

## Menizip

All-Inside Meniscal Repair System

CAT NO.	Variant
4133	4 Loaded (2 Pairs)
4144	2 Loaded (1 Pair)

# Deepfix

Screw-In PEEK Suture Anchor

CAT NO.	VARIANT	SIZE
3561	Triple Loaded	5.5 mm
3562	Double Loaded	5.5 mm
3564	Double Loaded	4.0 mm

# CortiZip

Soft Tissue Repair PEEK Anchor

CAT NO.	VARIANT	SIZE
3563	–	5.5 mm

# Tenemorf RC

All-Suture Anchor

CAT NO.	VARIANT	SIZE
3567	Triple Loaded	2.8 mm
3568	Double Loaded	2.8 mm

# PeriZip

PEEK Suture Anchor

CAT NO.	VARIANT	SIZE
3565	Double Loaded	2.5 mm
3566	Single Loaded	2.5 mm

# Tenemorf

All-Suture Anchor

CAT NO.	VARIANT	SIZE
3569	Double Loaded	1.8 mm
3570	Single Loaded	1.8 mm

## Disposable Full Thread Arthroscopy Cannula

CAT NO.	SIZE
3860	Ø8 X 70 mm
3861	Ø6 X 70 mm









**bioradmedisys**™  
science for people

[bioradmedisys.com](http://bioradmedisys.com)

Biorad Medisys Pvt. Ltd.

**Corporate Office & Plant:**

Survey No.: 48/3 & 48/7, Pashan Sus Rd, Sus Village,  
Taluka Mulshi, Pune - 411021. MH, India. Phone: 020-35000200

Follow us on    

BRP/CAT/AR/ART/01/00; 07-23